



ACC # 726855  
DP-MS-78-85

A PULSE HEIGHT ANALYZER/OSCILLOSCOPE DISPLAY INTERFACE

SRL  
RECORD COPY

by

W. J. Woodward, Jr.  
Savannah River Laboratory  
Aiken, SC 29801

Proposed for publication in  
*Nuclear Instruments and Methods*

---

This paper was prepared in connection with work under Contract No. AT(07-2)-1 with the U. S. Department of Energy. By acceptance of this paper, the publisher and/or recipient acknowledges the U. S. Government's right to retain a nonexclusive, royalty-free license in and to any copyright covering this paper, along with the right to reproduce and to authorize others to reproduce all or part of the copyrighted paper.

This document was prepared in conjunction with work accomplished under Contract No. DE-AC09-76SR00001 with the U.S. Department of Energy.

## **DISCLAIMER**

This report was prepared as an account of work sponsored by an agency of the United States Government. Neither the United States Government nor any agency thereof, nor any of their employees, makes any warranty, express or implied, or assumes any legal liability or responsibility for the accuracy, completeness, or usefulness of any information, apparatus, product or process disclosed, or represents that its use would not infringe privately owned rights. Reference herein to any specific commercial product, process or service by trade name, trademark, manufacturer, or otherwise does not necessarily constitute or imply its endorsement, recommendation, or favoring by the United States Government or any agency thereof. The views and opinions of authors expressed herein do not necessarily state or reflect those of the United States Government or any agency thereof.

This report has been reproduced directly from the best available copy.

Available for sale to the public, in paper, from: U.S. Department of Commerce, National Technical Information Service, 5285 Port Royal Road, Springfield, VA 22161, phone: (800) 553-6847, fax: (703) 605-6900, email: [orders@ntis.fedworld.gov](mailto:orders@ntis.fedworld.gov) online ordering: <http://www.ntis.gov/ordering.htm>

Available electronically at <http://www.doe.gov/bridge>

Available for a processing fee to U.S. Department of Energy and its contractors, in paper, from: U.S. Department of Energy, Office of Scientific and Technical Information, P.O. Box 62, Oak Ridge, TN 37831-0062, phone: (865 ) 576-8401, fax: (865) 576-5728, email: [reports@adonis.osti.gov](mailto:reports@adonis.osti.gov)

## ABSTRACT

A portable interfacing instrument was developed that accepts digital information from a pulse height analyzer and displays it as an energy spectrum on an ordinary laboratory oscilloscope.

## A PULSE HEIGHT ANALYZER/OSCILLOSCOPE DISPLAY INTERFACE

### INTRODUCTION

A 4096-channel Pulse Height Analyzer (PHA) and a Microprocessor-based Controller were developed at the Savannah River Laboratory.<sup>1,2</sup> An energy spectrum stored in a memory array in the PHA may be read out as 4096 binary words, each 20 bits long. All or part of this spectrum would best be displayed in graphic form as discrete energy levels on the horizontal axis versus total number of events at each level on the vertical axis. Ordinary laboratory oscilloscopes are available for this type display, therefore an instrument was developed to convert the digital PHA data to the required oscilloscope input signals.

Initial design requirements were:

1. The device must be portable, and therefore light weight, a convenient size, and battery-operated. (rechargeable batteries would be preferred.)
2. Vertical display must be linear, the analog equivalent of 10 successive binary bits, and range selectable from the 10 least significant through the 10 most significant bits in 6 steps.

3. Horizontal display must be linear, include 4096 channels (or less, as desired), and be a function of the timed sweep of the oscilloscope.

#### COMPONENT SELECTION

The design requirements dictated the choice of most of the system components. Some are listed below:

1. COS/MOS (Complimentary-symmetry, metal-oxide semiconductor) logic for its low power consumption.
2. Sealed lead-acid, gelled-electrolyte batteries for their capacity, charge-discharge characteristics, and convenient voltage ratings.
3. Wire-wrapped logic boards for their ease of plugging in dual in-line packaged (DIP) chips, headers, and extra-board ribbon cable connectors.
4. An aluminum attaché-sized case for its ease of handling, light weight, and durability.

The choice of COS/MOS logic influenced the selection of the Digital-to-Analog Converter (DAC). In addition to COS/MOS compatibility, the DAC must have good resolution and linearity plus a settling time fast enough to convert 4096 data words in a single oscilloscope sweep where the sweep repeats rapidly enough to provide a flicker-free display. Analog Devices' model AD 7520 meets

these requirements, will operate from a single 5 to 15 v power supply, and dissipates only 20 MW of power. Output is derived from a separate reference voltage divided via a switched R-2R ladder network. The DAC is a 10-bit unit, of which the least significant bit represents less than 0.1% of full scale. Therefore, readability on an oscilloscope screen will not suffer even when the 10 bits converted are the most significant of the 20 available for display.

#### PULSE HEIGHT ANALYZER MEMORY ACCESS

The PHA memory is arranged so the 4096 locations, or addresses can be accessed serially when using the "DISPLAY" operating mode. A train of 4096 evenly-spaced pulses may be fed to an address scaler to generate 4096 12-bit addresses. These are applied sequentially to the 20 memory chips ( $4096 \times 1$  bit each), each of which goes to a readout buffer. The square-wave pulse train causes the address scaler to advance on each rising edge. New data is latched into the buffers approximately 1  $\mu$ sec after each falling edge. Thus, 4096 words are presented serially at the buffer outputs (20 lines), each word a 20-bit binary number.

#### HORIZONTAL DISPLAY

A 250 KHZ clock frequency was chosen to operate the readout logic because of the 1  $\mu$ sec data switching time of the memory, the 500 nsec settling time of the DAC, and the sweep time

of approximately 16 msec for 4096 channels. This requires a sweep repetition rate of a little less than 60 Hz, which yields a flicker-free oscilloscope trace. The 250 KHz pulse train is easily obtained by counting down from the 4 MHz clock frequency in the PHA, and is presented to the interface as a square wave between high and low logic levels - a continuous train of pulses 2  $\mu$ sec wide and 2  $\mu$ sec apart.

The portion of the instrument panel (Fig. 1) marked "SCOPE DISPLAY" contains the mode switch and two 4-decade digital thumb-wheel switches (Fig. 2). On the panel, the mode switch is marked "HORIZONTAL." When in the "SET UP" position, the mode switch selects gating and timing logic to enable a 4-digit counter ("B") to count the 250 KHz clock pulses until 4096 have been registered, at which time the counter is disabled. 4096 pulses are gated out (4KCLK) to the PHA and the scope vertical input simultaneously. The scope sweep is triggered at the beginning of 4KCLK. The operator may use this display to adjust the scope's sweep controls so that the 4096 pulses just fill the desired portion of the screen. A timer in the interface is started when the scope is triggered, and the time may be adjusted to allow for full sweep and retrace, after which the process is repeated. The timer adjustment is marked "STABILITY," and is used to eliminate display flicker. The PHA address scaler is reset after the 4096<sup>th</sup> pulse via the clear display (CLRDPY) line; counter "B" is reset at the same time.

The same gating/timing action takes place when the mode switch is in the "4096 CHANNELS" position, except that the scope vertical input is connected to the analog output of the DAC instead of 4KCLK.

Positions 3 and 4 of the mode switch bring the digital thumb-wheel switches into play. Panel markings read "INTENSIFY (or EXPAND) XXXX CHANNELS STARTING WITH CHANNEL NO. XXXX," where XXXX is a 4-digit decimal number set up on the thumb-wheels by the operator. The "INTENSIFY" position allows display of all 4096 channels (sweep triggered as before), but the demultiplexer outputs are inhibited for portions of alternate sweeps. On odd-numbered sweeps analog output is zero from channel "0" through the number set up for counter "B", normal for XXXX (B) through  $[XXXX (B) + XXXX (A)]$ , then zero again through channel 4095. On even-numbered sweeps analog output is normal from 0 through 4095. The display intensity may be adjusted so that a difference in intensity and flicker is noticeable between the region of interest and the remainder of the sweep. This allows the operator to locate a peak in the spectrum by using the thumbwheel switches. The sweep may be calibrated in terms of energy for further peak identification.

In the "EXPAND" mode the sweep is triggered at the beginning of the region of interest (XXXX for counter "B"), and the sweep speed may be increased to expand that portion of the spectrum as desired.



## VERTICAL DISPLAY

The range switch shown in Figure 2 is marked "VERTICAL" (see Figure 1) and has 8 positions. This switch sends binary-coded control information to 10 one-of-eight demultiplexers (DEMUX) as shown in Figure 3. The first two inputs to the 10 DEMUX chips are connected to logic low and high, respectively. The first two range positions ("CAL: ZERO, FULL SCALE") connect the 10 DAC inputs to these for calibration; i.e., all bits low yield zero analog output and all bits high yield full scale output. The scope vertical gain may be adjusted as desired using these extremes for reference.

The remaining switch positions connect the DAC (via the DEMUXs) to 10 data input bits at a time ( $DI_0$  thru  $DI_9$ ,  $DI_{12}$  thru  $DI_{11}$ , ---- $DI_{10}$  thru  $DI_{19}$ ) in 6 steps. Switch positions are marked "COUNTS FULL SCALE: 1K, 4K, ---- $10^6$ ." It should be noted that counts in excess of full scale may appear on screen as a split image. Adjusting the range switch will correct this if the total does not exceed the maximum range.

## BATTERIES

An instrument was also designed for recharging the interface unit's batteries (Figure 4). Charging is accomplished by feeding a constant current of 150 ma to the batteries until their series voltage reaches 14.4 v, then switching to constant 25 mA to the batteries and an indicator diode in parallel. As battery charge increases the float voltage of 13.2 is reached and maintained indefinitely.

The power switch (Fig. 1) has 4 positions: "OFF," where batteries and external source are disconnected; "INT. BATT.," where only the batteries are connected to the load; "EXT. CHG.," where batteries and load are connected to the external charging source; and "CHG. ONLY," where the load is disconnected but the battery is connected to the charging source. A small light-emitting diode indicates when the external charger is supplying current, and a meter indicates battery terminal voltage in all but the "OFF" switch positions. An auxiliary switch disconnects the batteries when the instrument case is closed to prevent inadvertent discharge.

The two 6-v batteries are series connected, and are easily accessible behind a removable cover on the panel.

#### PERFORMANCE OF THE INTERFACE DEVICE

The prototype PHA/oscilloscope interface has been operated at SRL in laboratory and field service for approximately one year. Although operation has been satisfactory, it has been noted that switching noise in the DAC's analog output creates a background glow on the oscilloscope screen at normal brightness settings. This glow appears as a band about the viewed spectrum, the width of which depends on scope vertical gain and interface range settings. The background band does not diminish effectiveness or readability of the display; however, future models will be designed to eliminate this noise.

Also, when the displayed spectrum represents near zero on the vertical axis, the "INTENSITY" mode of operation loses its effectiveness because alternate sweeps cover essentially the same trace. This does not present a problem because that mode of operation is not normally used except when it is desired to locate one region of interest amid other data-containing channels.

The instrument has been used with Phillips, Telequipment, and Tektronix oscilloscopes. Its operation has been trouble-free.

## REFERENCES

1. J. S. Byrd and M. H. Goosey, "A Pulse Height Analyzer with Automatic Deadtime Correction," DP-MS-78-19. (Submitted for publication in *Nuclear Instruments and Methods*.)
2. J. S. Byrd, "A Microprocessor-based Controller for Nuclear Instruments." *IEEE Transactions on Nuclear Science*, Oct. 1978.



FIGURE 1. Instrument Front Panel

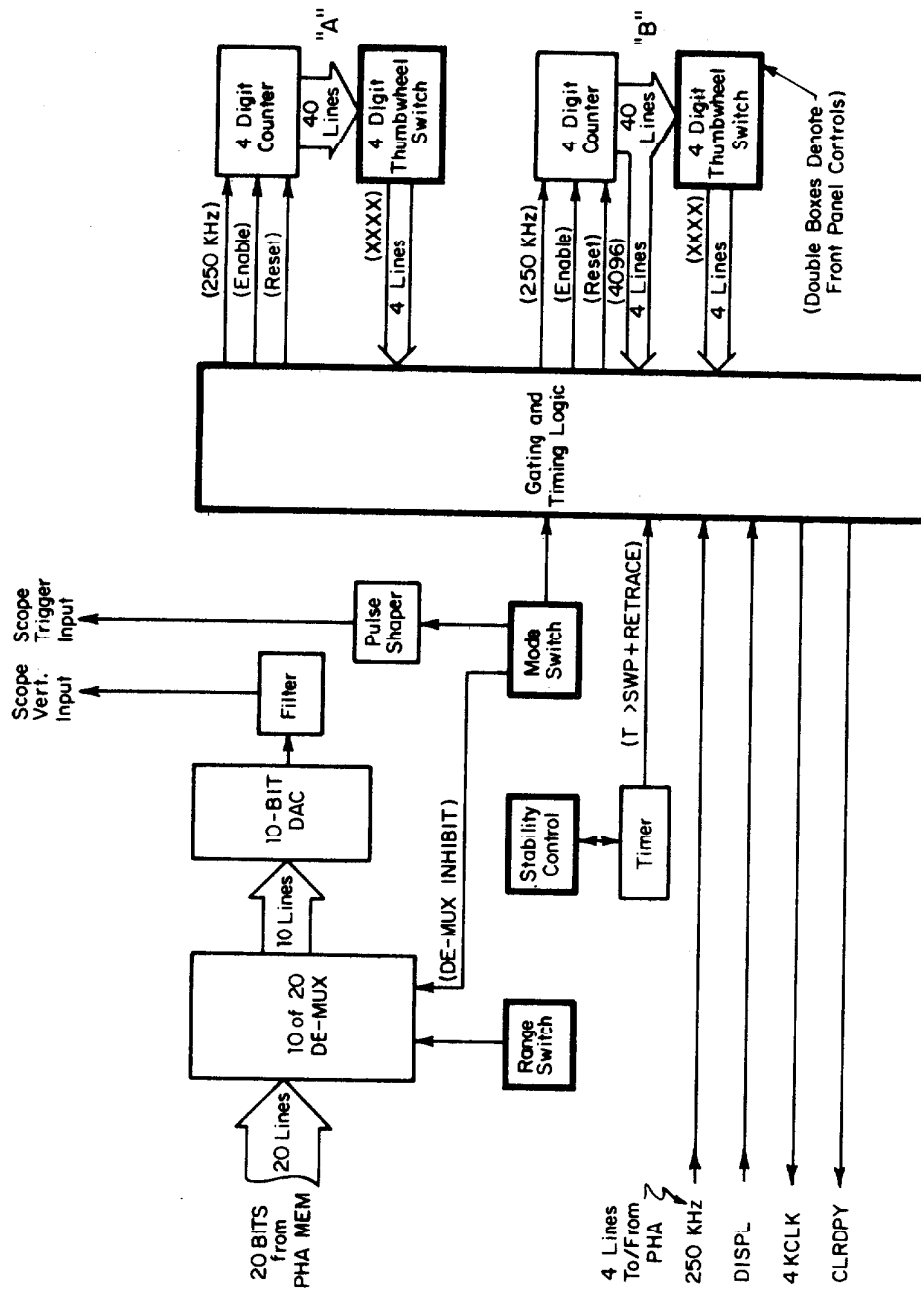


FIGURE 2. PHA/Oscilloscope Interface, Simplified Block Diagram

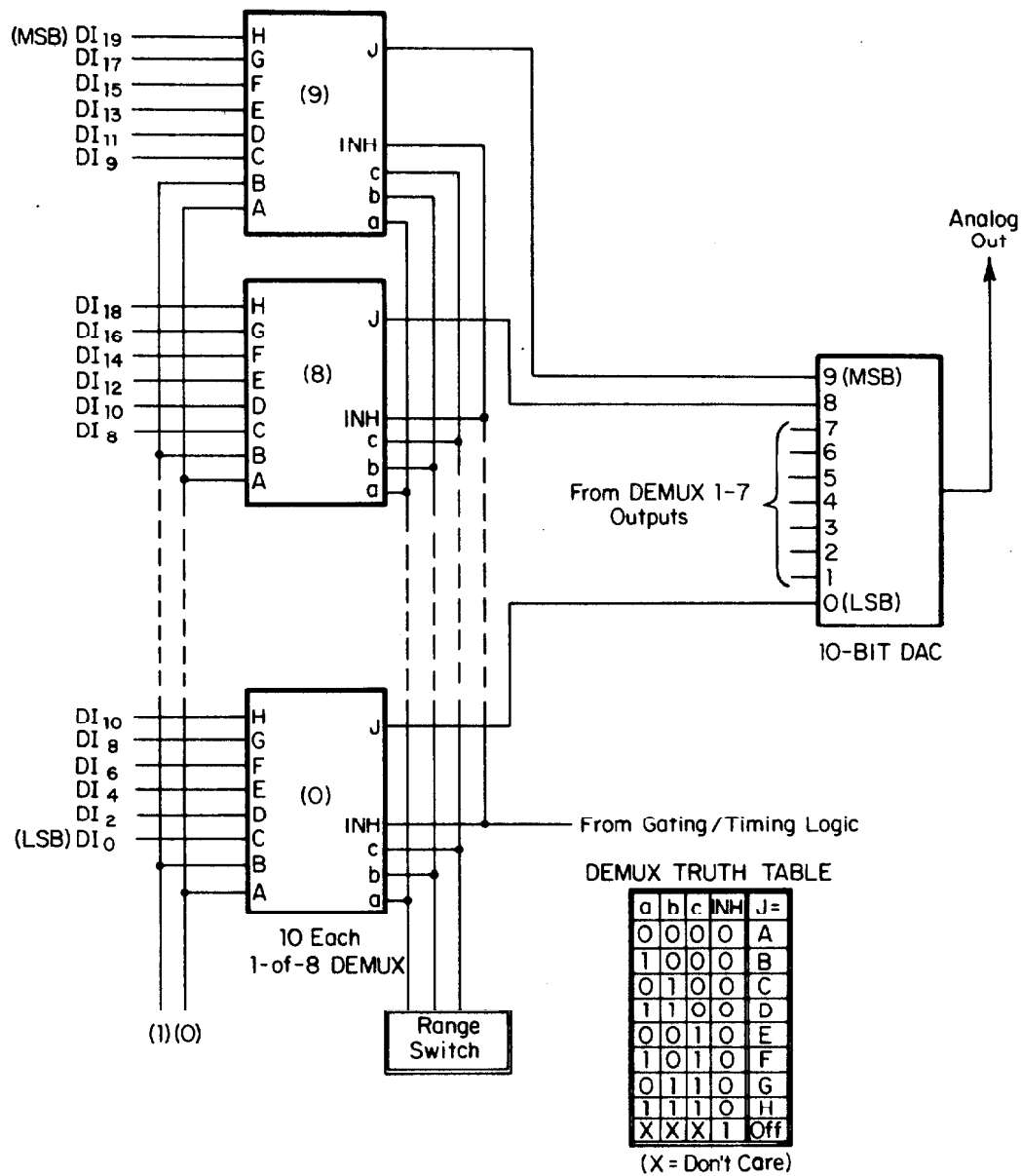


FIGURE 3. Demultiplexer-DAC, Simplified Diagram

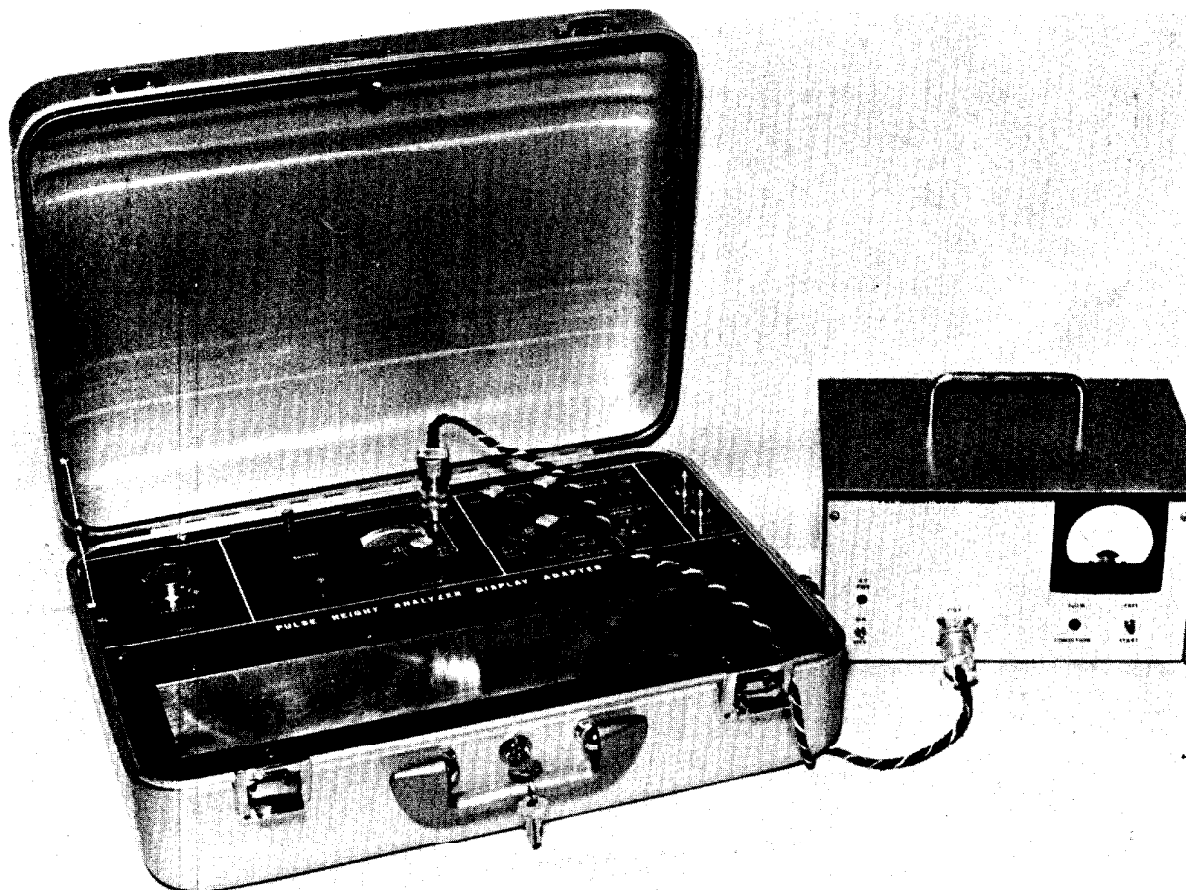


FIGURE 4. Instrument with Battery Charger Connected